

# CG-900 Cryo-Blade™

Cryogenic grinder for high-throughput QuEChERS, food safety and pesticide extraction applications



Breakthrough, automated processing from bulk to test portions in 15 minutes or less

- **High-throughput** – process up to 1 kg of sample per run
- **Efficient** – complete a method in 15 minutes or less
- **Versatile** – capability to use liquid nitrogen or dry ice to process your sample
- **Safety** – automated dispensing of liquid nitrogen eliminates the handling and manipulation of liquid nitrogen (patent pending)

The Cryo-Blade is an ultra high-capacity cryogenic grinder for comminution and processing capability up to 1 kg. The homogeneity of the sample is very efficient. Only 1 g is needed as a representative test sample from a 1 kg ground batch. The high-throughput capability saves valuable time for your laboratory and meets Food and Drug Administration (FDA) regulations. The Cryo-Blade automatically dispenses liquid nitrogen to allow the sample to reach cryogenic temperatures and then activates a grinding phase to pulverize samples, resulting in repeatable and reproducible sample preparation. No liquid nitrogen handling required. When liquid nitrogen is not desired or available, the unit can be used with dry ice to achieve desired results.

Automatic and programmed modes are available. In automatic mode, parameters are set based on the volume of the sample. Programmed mode allows users to save up to 10 methods. The intuitive interface includes 12 languages.

## FEATURES

- Intuitive color LCD screen interface
- Patented automatic dispensing of liquid nitrogen
- Complete a method in 15 minutes or less
- High-throughput processing up to 1 kg

## APPLICATIONS INCLUDE

- Analysis of organic pollutants in food products (pesticide residues)
- QuEChERS
- The unit accommodates an extensive range of both dry and wet samples from watermelon to the toughest of samples such as muscle tissue

## SAMPLE TYPES INCLUDE

- Fruits
- Vegetables
- Fungi
- Animal tissue
- Seeds
- Spices

## Accessories and What's Included

The insulated outer bowl, mechanical locking lever, and sensors provide the user with a safe, consistent, and reliable sample each use. The removable 5 L inner bowl can accommodate samples up to 1 kg and allows samples to be stored in the freezer, giving the user the option of using a second bowl to prepare a new batch. The outer bowl is insulated to allow safe touch at all times, although it is recommended to use protective gloves.


What's included: Cryo-Blade unit, insulated outer bowl, inner bowl, cryogenic processor lid with blade USB memory stick, and instruction manual.



CG-900 Cryo-Blade™  
Cryogenic Grinder



Cryogenic Process Inner  
Bowl, 5 L

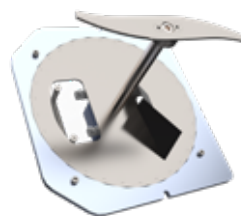

Cryogenic Processor Blade

## Before and After Samples



Watermelon



Broccoli



Fish



Ginger

For additional information or to request a quote or demo, scan the QR code or visit [cpsampleprep.com](https://cpsampleprep.com)



## Specifications

Specification	CG-900
Type of mill	Cryogenic grinder
Interface	Touchscreen, 12 languages
Motor	1 hp
Mixing speed	180 rpm
Max grinding speed	3600 rpm
Sample volume	1 kg
CE Approved	Yes
Voltage	115 VAC, 60 Hz or 230 VAC, 50 Hz
Power cord	3-prong grounded plug for 115 VAC or 2-prong European plug for 230 VAC
Dimensions (W x D x H)	45.0 x 63.0 x 70.0 cm (17.7 x 24.8 x 27.5 in)
Net weight	50 kg (110 lb), without grinding container
Warranty	1 year

## Ordering Information

Description	Item Number
CG-900 Cryo-Blade Cryogenic Grinder, 115 VAC, 60 Hz	CG-900-115
CG-900 Cryo-Blade Cryogenic Grinder, 230 VAC, 50 Hz	CG-900-230
Cryogenic Processor Inner Bowl, 5 L	41508
Cryogenic Processor Blade	41514

### Need a Demo?

**Virtual sample processing** – virtual demonstration of the equipment and sample processing

**Physical sample processing** – physical testing of your sample with shared results

**Physical equipment testing** – a demo unit can be sent to your facility to test before purchasing

Contact us today at [sampleprep@coleparmer.com](mailto:sampleprep@coleparmer.com) for details.

65 Liberty Street  
Metuchen, NJ 08840, USA  
Phone: 1.732.623.0465  
Email: [sampleprep@coleparmer.com](mailto:sampleprep@coleparmer.com)

**Cole-Parmer®**  
sampleprep  
[cpsampleprep.com](https://cpsampleprep.com)